

Product Data Sheet

DIAION™ SMT50

DIAION™ SMT50 is a mixed resin with strongly acidic cation exchange resin, and strongly basic anion exchange resin. It is used for non-regenerable mixed bed ion exchange applications for higher purity water.

Product

Grade Name	DIAION™ SMT50	
Type	Mixed	
Matrix	Styrene-DVB, Gel	
Functional Group	Sulfonic acid / Type I (trimethyl ammonium groups)	
Ionic Form	H ⁺ / OH ⁻	
Chemical Equivalent Ratio	1 / 1	

Specification

Component	Mixed resin	
Outlet Resistivity	MΩ·cm	18 min.
ΔTOC	ppb	5 max.

Typical Properties

Component	Mixed resin		
Shipping Density	g/L	710	
Component		Cation exchange resin	Anion exchange resin
Whole Bead Count	-	90 min.	-
Salt Splitting Capacity	meq/mL	1.7 min.	0.9 min.
Water Content	%	50 - 60	62 - 72
Particle Size Distribution on 1180 μm	%	5 max.	5 max.
Particle Size Distribution thr. 300 μm	%	1 max.	1 max.
Mean Particle Size	μm	700	720
Effective Size	mm	0.40 min.	0.40 min.
Uniformity Coefficient	-	1.6 max.	1.6 max.
Ionic Form Conversion (H ⁺)	eq%	99 min.	-
Ionic Form Conversion (OH ⁻)	eq%	-	90 min.
Ionic Form Conversion (Cl ⁻)	eq%	-	1 max.
Particle Density	g/mL	1.20	1.08
Total Swelling (Na ⁺ to H ⁺)	%	9	-
Total Swelling (Cl ⁻ to OH ⁻)	%	-	24

Recommended Operating Conditions

Maximum Operating Temperature	°C	60
Operating pH Range		0 - 14
Minimum Bed Depth	mm	800
Service Flow Rate	m/h	10 - 60

Hydraulic Characteristics

The approximate pressure drop at various temperatures and flow rates for each meter of bed depth of DIAION™ SMT50 resin in normal down flow operation is shown in the graphs below.

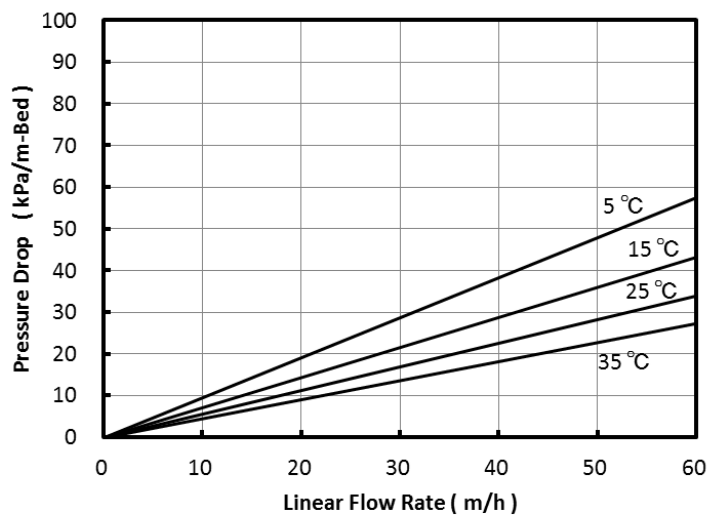


Fig. 1 Pressure Drop of SMT50

Rinse Performance

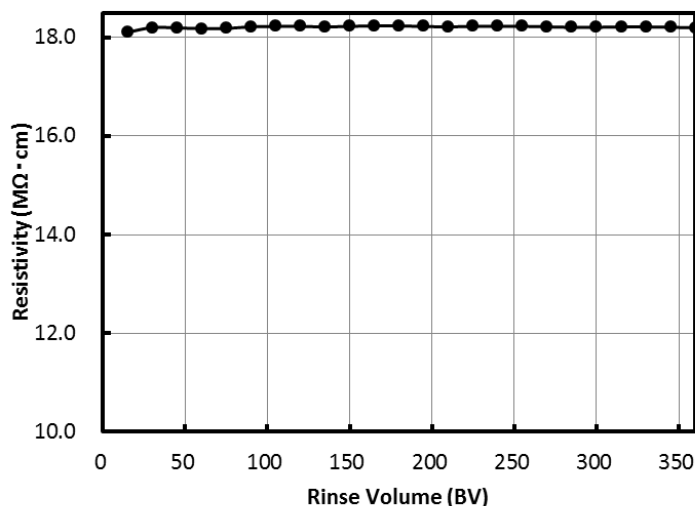


Fig. 2 Resistivity versus Rinse Volume for SMT50
Flow rate : SV 30 (15 L/hr), Resin volume : 500 mL-R

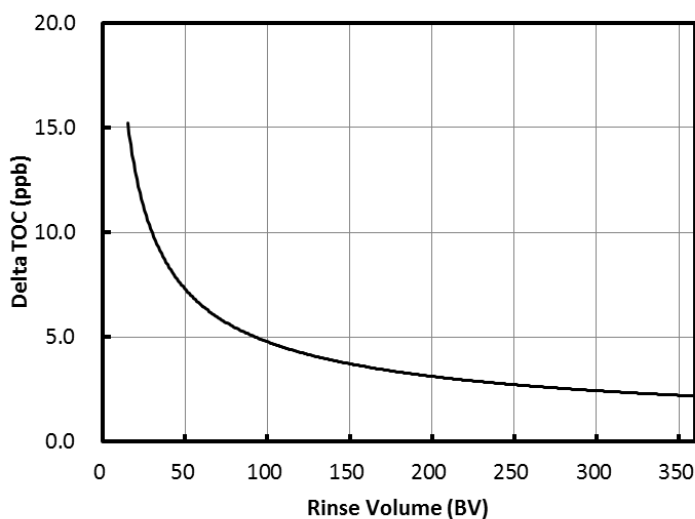


Fig. 3 Delta TOC versus Rinse Volume for SMT50
Flow rate : SV 30 (15 L/hr), Resin volume : 500 mL-R

Notice

This information are given in good faith but without warranty, and this also applies where proprietary rights of third parties are involved. The application, use and processing of our products are beyond our control and therefore your own responsibility.