

## Product Data Sheet

## DIAION™ WA30C

DIAION™ WA30C is a highly porous type weakly basic anion exchange resin. It has tertiary amine functionality and low-leachable properties compared with DIAION™ WA30. A wide range of applications, especially in a field of higher purity water treatment application, is recommended.

## Product

Grade Name	DIAION™ WA30C	
Type	Weak Base Anion	
Matrix	Styrene-DVB, Highly Porous	
Functional Group	Tertiary Amine	
Ionic Form	Free Base	

## Specification

Whole Bead Count	-	95 min.
Total Exchange Capacity	meq/mL	1.5 min.
Water Content	%	43 - 55
Particle Size Distribution on 1180 µm	%	5 max.
Particle Size Distribution thr. 300 µm	%	1 max.
Effective Size	mm	0.40 min.
Uniformity Coefficient	-	1.6 max.
COD	ppm	5 max.

## Typical Properties

Shipping Density	g/L	650
Mean Particle Size	µm	700
Particle Density	g/mL	1.05
Total Swelling (FB to Cl <sup>-</sup> )	%	23

## Recommended Operating Conditions

Maximum Operating Temperature	°C	100
Operating pH Range		0 - 9
Minimum Bed Depth	mm	800
Service Flow Rate	BV/h	10 - 40
Regenerant		NaOH
Regenerant Concentration	%	NaOH 1 - 4
Regenerant Level	% of ionic load	120
Regenerant Flow Rate	BV/h	2 - 6
Total Rinse Requirement	BV	5 - 10

1 BV(Bed Volume)=1 m<sup>3</sup>/m<sup>3</sup>-resin

## Hydraulic Characteristics

The approximate pressure drop at various temperatures and flow rates for each meter of bed depth of DIAION™ WA30C resin in normal down flow operation is shown in the graphs below.

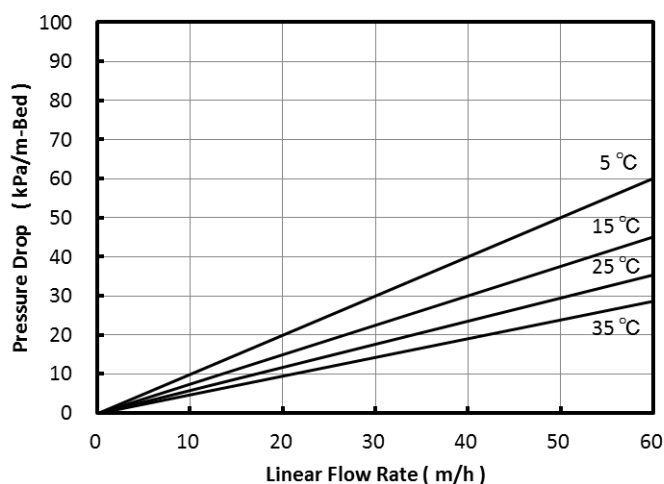


Fig. 1 Pressure Drop of WA30C

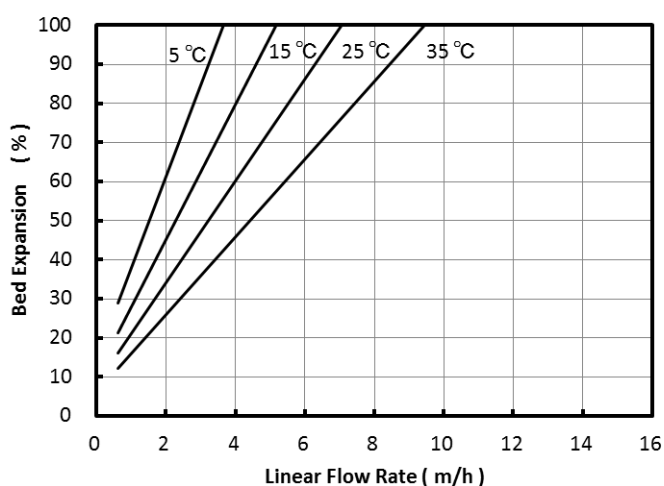


Fig. 2 Bed Expansion of WA30C

## Notice

DIAION™ is a registered trademark of Mitsubishi Chemical Corporation. The information contained herein is believed to be true and accurate, but all data, recommendations and suggestions are provided without guarantee, since the conditions of use are beyond our control and can affect the performance and properties of our products. The user is solely responsible for confirming that our product is suitable for the intended end use, and for compliance with all legal regulations and patents. Other than compliance with published Mitsubishi Chemical Corporation specifications agreed to pursuant to a signed writing during the warranty period, and except as required by law, MITSUBISHI CHEMICAL CORPORATION AND ITS AFFILIATES MAKE NO WARRANTIES, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OR ANY WARRANTY ARISING OUT OF A COURSE OF DEALING, CUSTOM OR USAGE OF TRADE. If a product is found to be defective during the warranty period, user's sole remedy and our sole obligation is, at our option, replacement of the affected product or refund of the purchase price. Except as required by law, we are not liable for any damage, harm or loss resulting from our product, whether direct, indirect, consequential, incidental or special, and irrespective of legal theory asserted, including strict liability, contract, warranty, or negligence.